

MagHoldIR Magnetic demountable liquid cell

Quick and easy sampling with disposable windows

Ideal for liquid samples with a volume as small as 1 μ L

The MagHoldIR is a sample holder using disposable windows suitable for micro-volume samples such as liquids, paints or adhesives. This holder can allow measurement of samples with a volume as small as 1 μ L. A sample is placed between the disposable KBr windows and secured by the magnetic cover. CaF₂ and BaF₂ windows are available as options.

Composition:

- Base plate (1 pc.)
- Cover plate (1 pc.)
- Micro KBr window (5 × 5 × 1 mm, 10 pcs.)
- Micro pipette (1 pc.)
- Pipette tip (5 pcs.)


Ordering information:

- 2000-0100 MagHold IR Magnetic demountable liquid cell (KBr window)
- 2000-0104 MagCA7 CaF₂ window, 2pcs/set
- 2000-0105 MagBA7 BaF₂ window, 2pcs/set



MagHoldIR


The MagHoldIR is a demountable sample holder for transmission measurement of micro samples such as liquids, paints or adhesive samples. The MagHoldIR is a magnetic sample holder that can provide measurement capability of samples with a volume as small as 1 μ L.

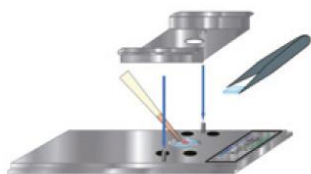

Samples are contained between KBr windows using only a magnetic cover plate.



Disposable KBr windows eliminates any window cleaning requirements.

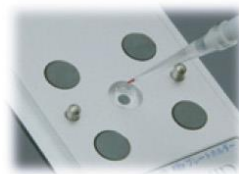

Optimized for micro measurements, small amounts of precious samples can be

Kit includes:	
MagHoldIR	× 1
Mini pipet	× 1
Pipet tips	× 5
Micro KBr windows	× 10

Replacement KBr : (2000-0047 KBr plate, 5x5x1mm, 100pcs./set)

JASCO France
17 Rue des Cerisiers
91090 Lisses - France

Tél. 01 64 97 09 60 Fax: 01 85 09 90 76
email: jascofrance@jascofrance.fr