


SL-5 E-Z™ Flow Cell

Our Model SL-5 E-Z™ Flow Cell is available in both sealed and demountable liquid flow cell configurations. The SL-5 E-Z™ Flow is made from stainless steel and features a Delrin® optical alignment cavity which, in the demountable version, precisely positions either a sandwich cell or the windows, spacers and gaskets in proper alignment with the cell filling ports. The cell features removable and replaceable compression nut and ferrule assemblies that hold the tubing in place. The standard nut is Delrin® and the standard ferrule is Tefzel®. However, stainless steel nut and ferrule assemblies are available. The cell can be used with both 1/16" OD and 1/8" OD tubing, and tubing diameters can be changed simply by removing the ferrule and nut assembly and inserting an inexpensive replacement. Special order high pressure versions of this cell are also available.



0007-7325	SL-5 E-Z™ Flow Cell- Demountable with Compression Fittings	1/16 Inch1/8 Inch	N/A	N/A
0007A-7305W NaCl	SL-5 E-Z™ Flow Cell- NaCl, with Compression Fittings	1/16 Inch1/8 Inch	0.1 mm0.2 mm0.5 mm	
0007B-7308W KCl	SL-5 E-Z™ Flow Cell- KCl, with Compression Fittings	1/16 Inch1/8 Inch	0.1 mm0.2 mm0.5 mm	
0007C-7311W KBr	SL-5 E-Z™ Flow Cell- KBr, with Compression Fittings	1/16 Inch1/8 Inch	0.1 mm0.2 mm0.5 mm	
0007D-7314W CaF2	SL-5 E-Z™ Flow Cell- CaF2, with Compression Fittings	1/16 Inch1/8 Inch	0.1 mm0.2 mm0.5 mm	
0007E-7317W BaF2	SL-5 E-Z™ Flow Cell- BaF2, with Compression Fittings	1/16 Inch1/8 Inch	0.1 mm0.2 mm0.5 mm	
0007H-7320W Infrasil Quartz	SL-5 E-Z™ Flow Cell- Infrasil Quartz, with Compression Fittings	1/16 Inch1/8 Inch	0.1 mm0.2 mm0.5 mm	
0007J-7322W ZnSe	SL-5 E-Z™ Flow Cell- ZnSe, with Compression Fittings	1/16 Inch1/8 Inch	0.1 mm0.2 mm0.5 mm	

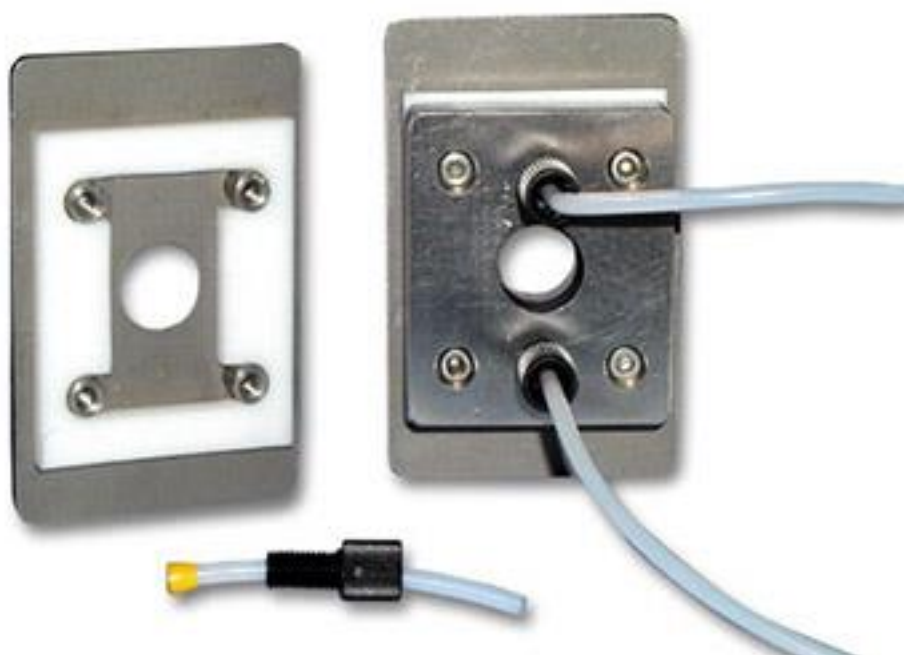

JASCO France
17 Rue des Cerisiers
91090 Lisses - France

Tél. 01 64 97 09 60 Fax: 01 85 09 90 76
email: jascofrance@jascofrance.fr